

# HiEpiXL™ Mammary Epithelial Cell Expansion Medium, Serum free

**Product Code: AL529**

## Product description:

HiEpiXL™ Mammary Epithelial Cell Expansion Medium is designed for *in vitro* cultivation and expansion of human mammary epithelial cells (HMEC) under serum free conditions. It contains basal medium (Part A), mammary epithelial cell growth supplement (Part B) and bovine pituitary extract (Part C). Part A consists of inorganic, organic salts, amino acids, vitamins and sodium bicarbonate. Part B consists of growth factors and nutrients necessary for growth of mammary epithelial cells. Part C is a bovine pituitary extract. This medium and supplement is devoid of antibiotics and antimycotics.

## Products Required But Not Supplied

1. Media Supplements	Code
Antibiotic-Antimycotic Solution 100X [or] Gentamicin-Amphotericin B solution 1000X	A002 A031
2. Reagents for Sub-culture	Code
Dulbecco's Phosphate Buffered Saline (DPBS)	TL1006
Trypsin/EDTA Solution 1X	TCL007
Trypan Blue 0.5% solution	TCL005
Trypsin inhibitor from soyabean	TCL068
3. Reagent for Coating Culture Vessel	Code
1% Collagen Solution in DPBS	TCL127

## Directions:

1. Thaw mammary epithelial growth supplement (Part B) and bovine pituitary extract (Part C) overnight at 2-8°C. *Note: Haziness in Part C after thawing is normal. This will not affect the performance of the medium.*

2. Disinfect the external surface of the bottles of Part A, Part B and Part C by spraying with isopropyl alcohol before placing in a biosafety hood.
3. Transfer the entire content of Part B and Part C to basal medium (Part A) under aseptic condition.  
*Note: If desired, 5ml of antibiotic - antimycotic solution (A002) or 0.5ml of gentamicin - amphotericin B solution (A031) can be added to 500ml of complete medium.*
4. Tightly cap the bottle and swirl gently to ensure proper mixing.  
*Note: Do not mix vigorously. Doing so will cause formation of foam.*
5. Store the complete medium at 2 - 8°C until use.

## Quality control:

### Appearance

Part A: Light pink coloured clear solution

Part B: Clear solution

Part C: Amber coloured solution

### pH of Part A

7.00-7.60

### Osmolality in mOsm/Kg H<sub>2</sub>O of Part A

320.00-360.00

### Sterility

No bacterial or fungal growth is observed after 14 days of incubation, as per USP specification.

### Cultural Response

The medium is tested for optimal cell growth and proliferation of human mammary epithelial cells.

### Storage and shelf life:

Store basal medium (Part A) at 2-8°C away from bright light. Store mammary epithelial cell growth supplement (Part B) and bovine pituitary extract (Part C) at -20 °C. Use before expiry date given on the product label. Shelf life of the complete medium is 4 weeks at 2-8°C.


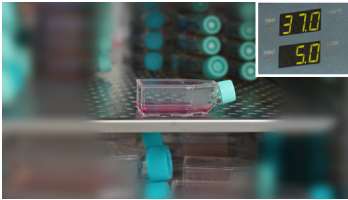

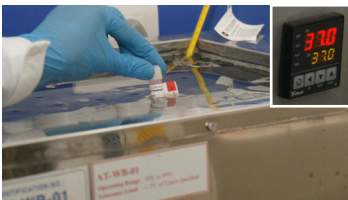


**Note:** Freezing of the basal medium and complete medium is not recommended. Avoid repeated freezing and thawing of the growth supplement.

Table 1 : Collagen Coating of Culture Vessel			
		Key Points to Remember	Time Required (approx.)
		<b>For uniform coating , make sure that the incubator is properly levelled</b>	
Aseptically add 1% collagen solution (TCL127)		<b>Refer Table 2 for recommended volumes of collagen solution</b>	1 min
Incubate for 2 hrs at 37°C incubator			2 hrs
Aspirate collagen solution with the help of pipette			
<b>COLLAGEN COATED CULTURE VESSEL IS READY FOR USE</b>			
If the vessel is not to be used immediately, store at 2-8°C upto one week.		<b>Flasks should be kept with caps tightly closed and plates should be sealed with a parafilm during storage</b>	

Table 2 : Recommended Volumes of Collagen Solution for Different Culture Vessels	
Culture Vessel	Volume Per Well
96-well plate	75 µl
48-well plate	150 µl
24-well plate	300 µl
12-well plate	500 µl
6-well plate	1 ml
T-25 Flask	5 ml
T-75 Flask	10 ml

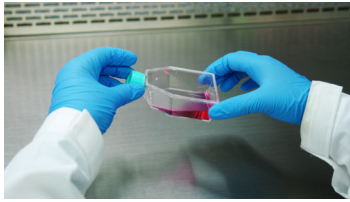
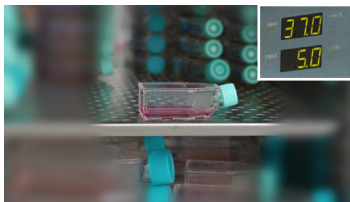
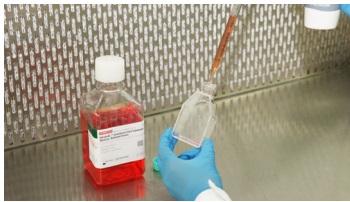
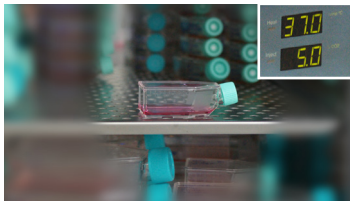
**Table 3 : Protocol for Thawing**

- Cryopreserved cells are supplied in liquid nitrogen dry vapour shipper (-150°C to -130°C).
- Upon receipt, immediately transfer the vial to the vapor phase of liquid nitrogen tank.
- Store it in the tank until further use. Cells must be processed at least in a BSL II hood.

		Key Points to Remember	Time Required (approx.)
<b>1. Preparation of Culture Vessel</b>			
a. Add 5ml of complete medium to a T-25 flask		<b>Preparation of complete medium</b> AL529 (Part A-500 ml) + (Part B-1.2 ml) + (Part C-2 ml) + A002(5ml)	60 secs
b. Place the flask at 37°C to equilibrate the medium			30 mins
<b>2. Thawing Procedure</b>		<b>Make sure water bath is set at 37°C before starting the thawing procedure</b>	
a. Remove cryovial from the liquid nitrogen tank/ shipper wearing appropriate protective gear		Thawing should be <b>AS FAST AS POSSIBLE</b> to minimize cell damage	
b. Immediately thaw the vial partially by holding in a water bath at 37°C		<b>DO NOT hold the vial</b> in water bath for more than 90-120 secs <b>AVOID</b> getting water upto the cap of the vial	90-120 secs
c. Disinfect the vial by swabbing thoroughly with 70% isopropyl alcohol			10 secs
d. Add the cell suspension drop by drop to the T-25 flask containing the pre-warmed complete medium. Keep the flask swirling while adding the cell suspension		<b>Dropwise addition</b> is required to prevent the cells from stress induced by exothermic reaction	30-60 secs


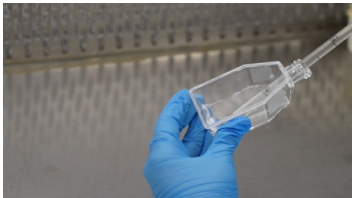

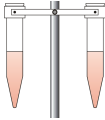

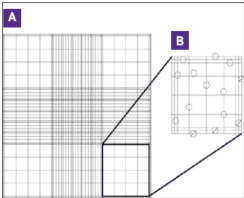
**Table 3 : Protocol for Thawing**

- Cryopreserved cells are supplied in liquid nitrogen dry vapour shipper (-150°C to -130°C).
- Upon receipt, immediately transfer the vial to the vapor phase of liquid nitrogen tank.
- Store it in the tank until further use. Cells must be processed at least in a BSL II hood.

		Key Points to Remember	Time Required (approx.)
e. Cap the flask and shake gently to ensure proper mixing and uniform distribution of cells in the medium			10 secs
<b>3. Incubation</b>			
a. Incubate the cells at 37°C and 5% CO <sub>2</sub>		Check for cell attachment in 2-3 hrs	2-3 hrs
b. If more than 70-80% cells are attached, replace the medium with fresh medium		<b>Medium change after 2-3 hours is mandatory to remove traces of DMSO</b> If cells have not attached, centrifuge the cell suspension at 1000 rpm for 7-8 mins, resuspend and seed in fresh medium	60-120 secs 7-8 mins
c. Incubate the cells at 37°C and 5% CO <sub>2</sub>			3-5 days
<b>YOUR CELLS ARE READY TO SUB-CULTURE</b>			
<b>4. Maintenance</b>			
a. Monitor the cells every day b. Change the medium c. Sub-culture once cells reach 70 - 80% confluence		Use the recommended freezing medium for cryopreservation of cells <b>Upto 50% Confluency:</b> <b>Change the medium on alternate day</b> <b>After 50% Confluency:</b> <b>Change the medium everyday</b>	

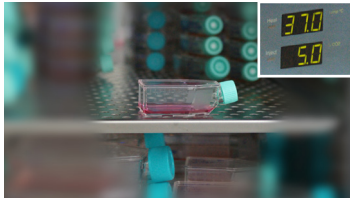
**Table 4 : Sub-culture**

- HMEC can be sub-cultured at a seeding density of 5000-10,000 cells/cm<sup>2</sup>.
- Sub-culturing ratios can vary from 1:2 - 1:5
- A confluent T-25 flask of HMEC yields 1.0 x 10<sup>6</sup> cells

		Key Points to Remember	Time Required (approx.)
a. Aspirate entire medium and discard DO NOT disturb the monolayer			60 secs
b. Wash the cells with 2-3 ml DPBS to remove residual medium c. Aspirate off the DPBS and discard		Prior to use, make sure that <b>Trypsin-EDTA solution is equilibrated to room temperature.</b>	60 secs
d. Add 2 ml pre-warmed Trypsin-EDTA (TCL128) solution e. Incubate the flask at 37°C for 1 min		Gently rock the flask to ensure complete coverage of the Trypsin-EDTA solution over the cells  <b>Exposing the cells to Trypsin-EDTA for longer time leads to loss of cell viability</b>	1 min
f. Check for rounding of the cells and keep tapping the flask gently for 1 min		This will allow <b>complete dissociation of cells</b>	1 min
g. To neutralize action of trypsin add 4 ml of TCL068 h. Transfer the entire contents of flask into 15ml centrifuge tube. i. Centrifuge the cells at 1000 rpm for 10 min		Very <b>small pellet</b> will be observed	15 secs 10 min
j. Carefully discard the supernatant by aspiration k. Add 2 ml of AL529 l. Pipette gently to get a homogenous mixture		<b>DO NOT POUR</b>  <b>Vigorous pipetting will stress the cells</b>	60 secs
m. Count cells using hemocytometer n. Seed at recommended seeding density in a new flask containing fresh complete medium Refer to <b>Table 5</b>		<b>DO NOT refrigerate cells after splitting</b> <b>Seed immediately</b>	10-15 mins

**Table 4 : Sub-culture**

- HMEC can be sub-cultured at a seeding density of 5000-10,000 cells/cm<sup>2</sup>.
- Sub-culturing ratios can vary from 1:2 - 1:5
- A confluent T-25 flask of HMEC yields 1.0 x 10<sup>6</sup> cells

		Key Points to Remember	Time Required (approx.)
o. Incubate in a humidified incubator at 37°C and 5% CO <sub>2</sub>			48 hrs
<b>Maintenance</b>			
a. Monitor the cells every day b. Change the medium c. Sub-culture once cells reach 70 - 80% confluence		<b>Upto 50% Confluency:</b> <b>Change the medium on alternate day</b> <b>After 50% Confluency:</b> <b>Change the medium everyday</b>	

**Table 5 : Seeding Density**

Flask	Recommended Seeding Density	No. of Cells Per Flask	Volume of Medium (ml)
T-25	5000 cells/cm <sup>2</sup>	0.125 x 10 <sup>6</sup>	5 - 7
	10,000 cells/cm <sup>2</sup>	0.25 x 10 <sup>6</sup>	5 - 7

These are recommended seeding densities from literature and our studies. Higher seeding densities do not cause any harm to the cells and reduce the required population doublings per passage. Lower seeding densities may cause cells to lose viability, detach during culture and in general take more population doublings to reach confluence.

## Related products:

Product name	Code	Packing
HiFi™ Human Mammary Epithelial Cells	CL014-0.5	0.5 million cells/vial
	CL014-2x0.5	2x0.5 million cells/vial
	CL014-T25	1 T25cm <sup>2</sup> flask
	CL014-2xT25	2 T25cm <sup>2</sup> flask
Trypsin 0.05% Solution 1X	TCL132-5X100ML	5 x 100ml
	TCL132-2X500ML	2 x 500ml
	TCL132-6X500ML	6 x 500ml
Trypsin – EDTA Solution 1X	TCL128-5X100ML	5 x 100ml
	TCL128-2X500ML	2 x 500ml
	TCL128-6X500ML	6 x 500ml
Trypsin Inhibitor from soybean 1X; Liquid	TCL068-1X100ML	1 x 100ml
	TCL068-5X100ML	5 x 100ml
Dulbecco's Phosphate Buffered Saline	TL1006-5X100ML	5 x 100ml
	TL1006-2X500ML	2 x 500ml
	TL1006-6X500ML	6 x 500ml
	TL1006-18X500ML	18 x 500ml
	TL1006-1X1000ML	1 x 1000ml
Antibiotic Antimycotic solution 100X, Liquid	A002-5X20ML	5 x 20ml
	A002-5X50ML	5 x 50ml
	A002-5X100ML	5 x 100ml
Gentamycin Solution	A005-5X20ML	5 x 20ml
	A005-5X50ML	5 x 50ml
Gentamycin - Amphotericin B Solution 1000X	A031-5X20ML	5 x 20ml
	A031-5X50ML	5 x 50ml

Revision: 2 / 2017

### Disclaimer:

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ Publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal diagnostic or therapeutic use but for laboratory, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.